## MT 800

IB MT 800 E-V11 1420

| Description of this Thermometer         | EN |
|-----------------------------------------|----|
| ① ON/OFF button                         |    |
| ② Display                               |    |
| ③ Battery compartment cover (back side) |    |
| ④ Measuring sensor                      |    |

##### Important Safety Instructions

- Follow instructions for use. This document provides important product operation and safety information regarding this device. Please read this document thoroughly before using the device and keep for future reference.
- This device may be used only for measuring body temperature!
- The minimum measurement time until the beep is heard must be maintained without exception!
- Ensure that children do not use this device unsupervised; some parts are small enough to be swallowed.
- Do not use this device close to strong electromagnetic fields such as mobile telephones or radio installations. Keep a minimum distance of 3.3 m from such devices when using this device.
- Protect the device from impact and dropping!
- Avoid bending the tip more than 45°!
- Avoid ambient temperatures above 60 °C. NEVER boil this device!
- Use only the commercial disinfectants listed in the section «Cleaning and Disinfecting» to clean the device.
- We recommend this device is tested for accuracy every two years or after mechanical impact (e.g. being dropped). Please contact Micro-life-Service to arrange the test.

**WARNING:** The measurement result given by this device is not a diagnosis! Do not rely on the measurement result only.

Batteries and electronic devices must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

Read the instructions carefully before using this device.

Type BF applied part

##### Turning on the Thermometer

To turn on the thermometer, press the ON/OFF button ①; a short beep signals «thermometer ON». A display test is performed. All segments should be displayed.

The last measurement reading will be shown on the display ② automatically for 2 seconds with the «M» icon. Then at an ambient temperature of less than 32 °C, an «L» and a flashing «C» appear at the display field ②. The thermometer is now ready for use.

##### Function Test

Correct functioning of the thermometer is tested automatically each time it is turned on. If a malfunction is detected (measurement inaccuracy), this is indicated by «ERR» on the display, and a measurement becomes impossible. In this case, the thermometer must be replaced.

##### Using the Thermometer

Choose the preferred measuring method. When taking a measurement, the current temperature is continuously displayed and the «C» symbol flashes. If the beep is heard 10 times and the «C» is no longer flashing, the predictive end-temperature has been determined and the thermometer can be read now.

35 short beeps will sound when the temperature is higher than 37.5 °C in order to alert the patient that he/she may have fever. Reference: Oral temperature.

To achieve comparable results allow a 1 minute interval time between measurements.

To prolong the battery life, turn off the thermometer by briefly pressing the ON/OFF button ①. Otherwise the thermometer will automatically turn off after about 10 minutes.

##### Reviewing the last Measurement

Press ON/OFF button ① for 2 seconds, when turning on the thermometer. The last measurement will be displayed together with a «M» for memory. Press the ON/OFF button again to switch off the device.

##### Measuring methods / Normal body temperature

##### ► In the mouth (oral) / 35.5 - 37.5 °C

Position the thermometer in one of the two pockets under the tongue, to the left or right of the root of the tongue. The measuring sensor ④ must be in good contact with the tissue. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced by inhaled/ exhaled air.

##### Approx. measuring time: 10 seconds!

##### ► In the anus (rectal) / 36.6 - 38.0 °C

This is the most reliable measuring method, and is especially suitable for infants and small children. Carefully insert the measuring sensor ④ of the thermometer 2 to 3 cm into the anal aperture.

##### Approx. measuring time: 10 seconds!

##### ► In the armpit (axillary) / 34.7 - 37.3 °C

We recommend measuring temperature orally or rectally. If using axillary measurement, a minimum measurement time of 5-10 minutes is required regardless of the beep sound.

##### Cleaning and Disinfecting


Name: Isopropyl alcohol 70 %; immerse: max. 24 hours.

##### Battery Replacement

When the «▽» symbol (upside-down triangle) appears in the display, the battery is flat and needs replacing. To replace the battery remove the battery compartment cover ③ from the thermometer. Insert the new battery with the + at the top. Make sure you have a battery of the same type to hand. Batteries can be purchased at any electrical store.

| Technical Specifications                                                              |                                                                                                                                                                                      |
|---------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Type:</b>                                                                          | Predictive maximum thermometer                                                                                                                                                       |
| <b>Measurement range:</b>                                                             | 32.0 <span> </span> °C to 43.9 <span> </span> °C<br>Temp. < 32.0 <span> </span> °C: display «L» for low (too low)<br>Temp. > 43.9 <span> </span> °C: display «H» for high (too high) |
| <b>Measurement accuracy:</b>                                                          | ± 0.1 <span> </span> °C between 34 <span> </span> °C and 42 <span> </span> °C                                                                                                        |
| <b>Operating conditions:</b>                                                          | 10 - 40 <span> </span> °C; 15-95 <span> </span> % relative maximum humidity                                                                                                          |
| <b>Storage conditions:</b>                                                            | -25 - +60 <span> </span> °C; 15-95 <span> </span> % relative maximum humidity                                                                                                        |
| <b>Battery:</b>                                                                       | 1.5/1.55 V; LR41                                                                                                                                                                     |
| <b>Battery lifetime:</b>                                                              | approx. 4500 measurements (using a new battery)                                                                                                                                      |
| <b>IP Class:</b>                                                                      | IP67                                                                                                                                                                                 |
| <b>Reference to standards:</b>                                                        | EN 12470-3, clinical thermometers; ASTM E1112; IEC 60601-1; IEC 60601-1-2 (EMC); IEC 60601-1-11                                                                                      |
| <b>Expected service life:</b>                                                         | 5 years or 10000 measurements                                                                                                                                                        |
| This device complies with the requirements of the Medical Device Directive 93/42/EEC. |                                                                                                                                                                                      |
| Technical alterations reserved.                                                       |                                                                                                                                                                                      |
| <b>Guarantee</b>                                                                      |                                                                                                                                                                                      |

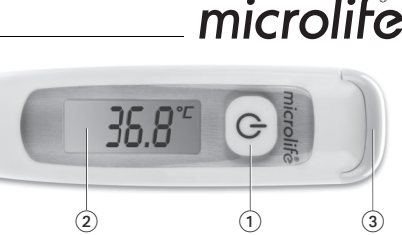We grant you a lifetime guarantee after the date of purchase. Any damage caused by improper handling shall not be covered by the guarantee. The battery and packaging are excluded from the guarantee. All other damage claims excluded. A guarantee claim must be submitted with the purchase receipt. Please pack your defective device well and send with official postage to your local Microlife-Service.

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