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Globalcare Medical Technology Co., Ltd
7th Building, 39 Middle Industrial Main Road,
European Industrial Zone, Xiaolan Town,
528415 Zhongshan City, Guangdong Province, PRC
Phone: +86 760 22589901
www.globalcare.com.hk/contact

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
EN → 1
RU → 6
KZ → 12

Microlife BC 100 Soft

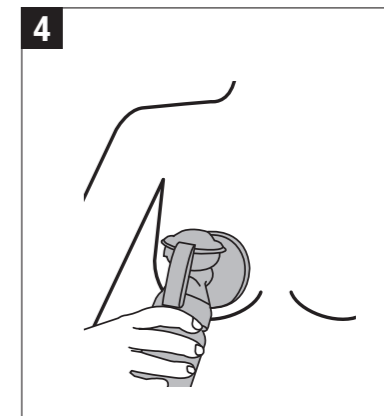- ① Manual pump attachment (2 levels)
- ② Silicone membrane
- ③ Silicone valve
- ④ Silicone cushion (standard size)
- ⑤ Bottle (180 ml)
- ⑥ Bottle holder
- ⑦ Teat
- ⑧ Bottle cap for teat
- ⑨ Screwed plug for bottle
- ⑩ Adapter for NUK (the device can also be used with this bottle if desired)
- ⑪ Adapter for AVENT (the device can also be used with this bottle if desired)

Read the instructions carefully before using this device.

NOTE:
Note on important information

WARNING
Warning instruction indicating a risk of injury or damage to health

IMPORTANT
Safety note indicating possible damage to the device/accessory

Unified sign of circulation of products on the market of the member states of the Eurasian Customs Union
Manufacturer

Serial number

Reference number

CE Marking of Conformity

Intended use:

This device must be used exclusively for pumping breast milk from nursing mothers. Do not use this device on animals. The device is only intended for domestic/private use, not for commercial use. Please read through these instructions carefully in order for you to understand all functions and safety information.

Dear Customer,

This device allows nursing mothers to express their breast milk easily and in the comfort of their own homes. If not enough breast milk is taken from mothers who have recently given birth, this can lead to the painful build-up of milk, which in the worst-case scenario could cause inflammation of the nipples (mastitis). This build-up of milk can be prevented through targeted expressing of breast milk using the breast pump. In addition, breast pumps help to alleviate cracked or sore nipples. The breast pump also allows you to store breast milk in the bottle, for example if you will be travelling or away from your baby for an extended period of time.

Babies who cannot be breast-fed directly for health reasons (due to a premature birth or cleft palate) can still be given breast milk by using the breast pump.

If you have any questions, problems or want to order spare parts please contact your local Microlife-Customer Service. Your dealer or pharmacy will be able to give you the address of the Microlife dealer in your country. Alternatively, visit the internet at www.microlife.com where you will find a wealth of invaluable information on our products.
Stay healthy – Microlife AG!

Table of contents

- 1. Important Safety Instructions
- 2. Using the device for the first time
- 3. Directions for Use
- 4. Malfunctions and Actions to take
- 5. Guarantee
- 6. Technical Specifications
- 7. Information about the sterility of the product
- 8. Disposal information

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1. Important Safety Instructions

Safety and protection

- Follow instructions for use. This document provides important product operation and safety information regarding this device. Please read this document thoroughly before using the device and keep for future reference.
- This device may only be used for the purposes described in these instructions. The manufacturer cannot be held liable for damage caused by incorrect application.
- Keep packaging material away from children to avoid danger of suffocation.
- Keep small parts out of reach of children.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- This device should only be used with original accessories as shown in these instructions.
- For hygiene reasons, the device may only be used on one person. Do not share the device with other nursing mothers.
- The device must only be used on the female breast.
- Never use the device if feeling drowsy or when asleep.
- Never use the device in the bath or the shower.
- Do not use this device if you think it is damaged or notice anything unusual.
- Never place the device in the microwave.
- Never heat up breast milk in the microwave as it may be heated unevenly and could burn your child when they drink it.
- This device must not be used by people with reduced physical, sensory or mental skills or a lack of experience or knowledge.
- Repairs must only be carried out by Customer Services or authorised retailers. Under no circumstances should you open or repair the device yourself, as faultless functionality can no longer be guaranteed thereafter. Failure to comply will result in voiding of the warranty.
- The device requires no calibration.
- No modification to the device is permitted.

WARNING

Child safety

- This device may only be used under adult supervision.
- Feeding from the teat continuously and for long periods can lead to tooth decay.
- Use the teat for feeding, not as a pacifier.
- Always check the temperature of the milk before feeding.
- Inspect the device before each use and pull the teat in all directions.
- Dispose of the teat at the first sign of damage or wear and tear.
- Do not expose the teat to direct sunlight or sources of heat and do not leave in disinfectant («sterilising solution») for longer than recommended as this may damage the teat.
- Keep the teat away from sweet substances as these may cause tooth decay.
- Do not use any abrasives or anti-bacterial cleaning agents.
- Sterilise all parts that come into contact with the breast milk for five minutes in boiling water.
- Clean the device before each use.
- Make sure that the temperature of the liquid in the bottle never exceeds 50 °C.

Ensure that children do not use this device unsupervised; some parts are small enough to be swallowed. Be aware of the risk of strangulation in case this device is supplied with cables or tubes.

IMPORTANT

When transporting breast milk in the bottle (e.g. when on the move or travelling), make sure that the bottle cap is always screwed on tightly.

2. Using the device for the first time

NOTE:

Check that the packaging of the device has not been tampered with and make sure that all components are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.

Assembling the device

Once the device has been removed from the packaging, check the plastic containers in particular for cracks.

1. Screw the pump attachment ① onto the bottle ⑤. Ensure the pump attachment is firmly secured so that no liquid can escape (see short instruction).
2. Attach the handle for the manual pump including the silicone membrane by pressing the silicone membrane firmly onto the pump attachment until it clicks into place (see short instructions).

Cleaning and sterilising the device

WARNING

Clean and sterilise the separate parts of the device that come into direct contact with the breast milk individually before first use and any further use.

To sterilise the device correctly, proceed as follows:

1. Unscrew the pump attachment ① from the bottle ⑤.
2. Remove the silicone membrane ② and silicone valve ③ from the pump attachment.
3. Sterilise the bottle and pump attachment (with silicone membrane and silicone valve) for five minutes in boiling water or using a steam steriliser. Do not use the device until all individual parts are completely dry. Do not use any antibacterial cleaning products, since these may damage the plastic on the device.

3. Directions for Use

Tips to ensure a better milk flow

- Choose a time and location to use the device where you can relax.
- Having your child nearby or looking at a picture of your child can have a positive effect on the milk flow.
- If feeding your child using one breast, use the device on the other breast at the same time.
- Warmth and relaxation encourage milk flow. Use the device after having a shower or bath.
- If your breasts are swollen or sore, place a warm towel on them for a few minutes. This will improve the milk flow.
- If the milk does not start to flow from your breasts immediately, try to relax and then try again a few minutes later.

- If the device still does not extract any milk from your breasts after five minutes, stop using it and try again later.
- If you experience any sharp pains when using the device, consult a doctor.

Operating the device

1. Before using the device, wash your hands and breasts thoroughly.
2. Sit with your upper body bent forwards slightly. You can support your back better by using a pillow.
3. Then place the breast shield with silicone cushion ④ gently on the breast. Make sure that the nipple is completely enclosed and there are no air pockets. To ensure proper function, the breast shield must be drawn firmly onto the breast.
4. You can vary the pump intensity on the manual pump attachment ① by selecting between the MIN/MAX markings that can be seen. We recommend starting in the MIN position.

NOTE:

See the pictures at the start of the instruction manual for more information.

Storing the breast milk

Expressed milk can be stored in the fridge or freezer for a certain period of time.

- When storing breast milk, ensure that the device and the bottles are disinfected before usage. Only store breast milk in disinfected containers.
- Store breast milk in the fridge or freezer immediately after expressing. Do not leave breast milk out at room temperature.
- If storing breast milk in the freezer, it is recommended that you mark the bottle with the date of expressing.

WARNING

Observe the storage periods in the following table:

	Room temperature	Refrigerator (approx. 3 - 5 °C)	Freezer (approx. -16 °C)
Recently expressed breast milk	Can be kept for 6 hours	Can be kept for 24 hours (Do not store in the fridge door, as it is not sufficiently cool there)	Can be kept for 3 months

	Room temperature	Refrigerator (approx. 3 - 5 °C)	Freezer (approx. -16 °C)
Defrosted breast milk (that had been frozen)	Use immediately	10 hours	Do not refreeze

Feeding with the breast milk

NOTE:

You can use the supplied bottle with the silicone teat when feeding with the breast milk. Alternatively, you can also express the breast milk with a NUK or Avent bottle before feeding (corresponding NUK / Avent adapter is included in the box).

- Prior to feeding, ensure the bottle has been disinfected before filling it with the breast milk that has been stored.
- To defrost breast milk, store it overnight in the fridge, so as to retain the nutrients. Once the breast milk is defrosted, it must be used within 24 hours. Breast milk should only be defrosted in warm water in especially urgent cases.
- If the breast milk smells unpleasant, it should be disposed of immediately.
- Never heat up breast milk in the microwave as it may be heated unevenly and could burn your child when they drink it.
- If heating up the breast milk, check its temperature before feeding it to your baby.

4. Malfunctions and Actions to take

Use following table for solutions if any problem occur.

Problem	Potential solutions
There is no suction from the device.	<ul style="list-style-type: none"> • Check all parts of the device are correctly assembled. • Check whether the pump attachment is correctly placed on the breast.
No breast milk is being pumped.	<ul style="list-style-type: none"> • Ensure that there is no milk remaining in the silicone valve and obstructing the milk flow.

Problem	Potential solutions
Pumping the breast milk is causing pain.	<ul style="list-style-type: none"> • Try to relax and then attempt it again a few minutes later. As you use the device more and more, you will find it easier to pump milk. • Set the suction level to MIN. • Consult your doctor.
The breast milk is leaking out underneath the silicone cushion.	<ul style="list-style-type: none"> • Remove the silicone cushion and re-attach it. Ensure the silicone cushion is correctly positioned and firmly secured to the pump attachment. • Lean your upper body forwards slightly.
One of the device components has become altered or damaged.	<ul style="list-style-type: none"> • Do not use any corrosive, gaseous cleaning agents or solvents. • In some cases, marked fluctuations in temperature can damage the plastic on the device. • If any of the parts becomes damaged, stop using the device.

5. Guarantee

This device is covered by a **2 year guarantee** from the date of purchase. During this guarantee period, at our discretion, Microlife will repair or replace the defective product free of charge.

Opening or altering the device invalidates the guarantee.

The following items are excluded from the guarantee:

- Transport costs and risks of transport.
- Damage caused by incorrect application or non-compliance with the instructions for use.
- Damage caused by accident or misuse.
- Packaging/storage material and instructions for use.
- Accessories and wearing parts: Pump attachment (with membrane/valve), tubes, cushion, bottle with holder, teat and cap, bottle adapter.

Should guarantee service be required, please contact the dealer from where the product was purchased, or your local Microlife service. You may contact your local Microlife service through our website:

www.microlife.com/support

Compensation is limited to the value of the product. The guarantee will be granted if the complete product is returned with the original invoice. Repair or replacement within guarantee does not prolong or renew the guarantee period. The legal claims and rights of consumers are not limited by this guarantee.

6. Technical Specifications

Type:	Manual Breast Pump BC 100 Soft
Dimensions:	200 x 160 x 66 mm
Weight:	135 g (without batteries)
Maximum suction pressure:	-0.32 bar
Minimum pumping pressure:	-0.08 bar
Operating conditions:	5 - 40 °C 30 - 85 % relative maximum humidity (non condensing) 700 - 1060 hPa ambient pressure
Storage conditions:	-20 - +60 °C 10 - 95 % relative maximum humidity (non condensing) 700 - 1060 hPa ambient pressure

- ▶ Child use and care articles - Drinking equipment according to EN 14350
- ▶ Device for expressing breast milk

The technical specifications may change without prior notice.

7. Information about the sterility of the product

The product is supplied non-sterile. Number of sterilization cycles: 360. Hand pump nozzle (2 levels), silicone membrane, silicone valve, silicone pad (standard size), bottle (180 ml), bottle nipple, bottle cap for nipple, bottle screw cap, adapter for NUK, adapter for AVENT steam sterilized before use for five minutes in boiling water or with a steam sterilizer.

8. Disposal information

Dispose of the product in accordance with SanPiN 2.1.3684-21 for class B waste (epidemiologically hazardous waste). After hardware methods of disinfection using physical methods and changing the appearance of waste, which excludes the possibility of their reuse, class B waste can be accumulated, temporarily stored, transported, destroyed and disposed of together with class A waste (epidemiologically safe waste, similar in composition to solid household waste). The packaging of class B decontaminated medical waste must be labeled indicating that the waste has been decontaminated.

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